


NEW

Any-Cem[®]

Self-adhesive Dual-cure Resin Cement for Zirconia and High-strength Indirect Restorations

Easy Application,
Strong Bond

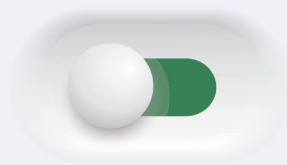

MDCLUS
MEDICLUS CO., LTD.

Any-Cem®

Fast and Durable Bonding for Zirconia

Any-Cem® is a self-adhesive resin cement that offers both self-cure and light-cure capabilities. Featuring a simplified application process and excellent mechanical bonding strength, it is recommended for the cementation of high-strength ceramic restorations, such as zirconia-based prosthetics. The formulation combines the advantages of both traditional cements and adhesive resin cements.

Streamlined Workflow with 10-MDP

Any-Cem® incorporates the specially developed functional monomer 10-MDP, which establishes stable chemical bonds with zirconia, various metal restorations, and tooth structures. This eliminates the need for additional surface pretreatment, enabling a faster and more efficient clinical workflow.


Outstanding bond strength

Unlike many competing products, Any-Cem® delivers high bond strength in both light-cure and self-cure modes. This ensures reliable adhesion even in areas with limited light access, such as metal-ceramic restorations. Its performance makes it also suitable for use with materials like **Titanium** and **Zirconia**.

Indications

- Zirconia: Crowns, Bridges, Endodontic posts
- All-ceramics: Crowns, 3-unit bridges
- Indirect composite resins: Endodontic posts
- Metal, Metal-ceramics: Inlays, Onlays, Crowns, Bridges, Endodontic posts
- Titanium: Implant abutment restorations

Step by Step Process

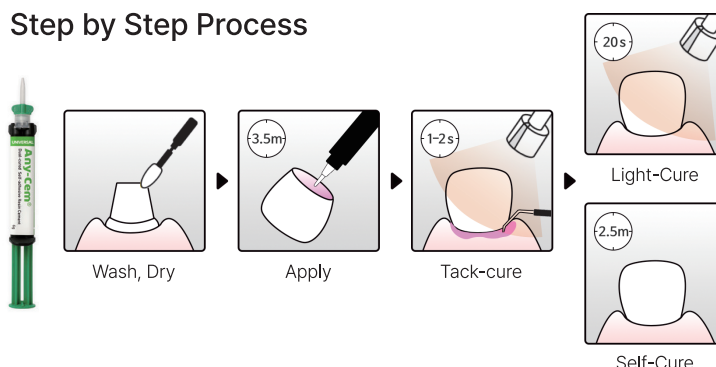

Benefits

- **Easy application**: No pretreatment or separate bonding required, thanks to 10-MDP
- **Efficient workflow**: No capsule mixing, specialized tools, or equipment needed
- **Versatile dispensing options**: Convenient Automix tips or economical Eco-tips
- **Effortless excess removal**
- **Convenient storage**: Room temperature (1–30 °C) stable

Delivery Forms



#203 Any-Cem® Universal

Any-Cem® 6g 1ea
Automix tips 10ea
Eco-tip 1ea
Mixing pad 1ea

#203-1 Any-Cem® White

Any-Cem® 6g 1ea
Automix tips 10ea
Eco-tip 1ea
Mixing pad 1ea

#203-3 Any-Cem® TL

Any-Cem® 6g 1ea
Automix tips 10ea
Eco-tip 1ea
Mixing pad 1ea

MEDICLUS

Headquarter | No.1210, 134, Gongdan-ro, Heungdeok-gu, Cheongju-si, Chungcheongbuk-do, 28576, Korea | +82 43 211 2877

Sales Office | No.206 Brown Stone, 464, Cheongpa-ro, Jung-gu, Seoul, 04510, Korea | +82 10 8385 2877

www.mediclus.co.kr | sales@mdclus.com

